

Methodology article

Open Access

Molecular beacon-based real-time PCR detection of primary isolates of *Salmonella* Typhimurium and *Salmonella* Enteritidis in environmental and clinical samples

Andreas V Hadjinicolaou^{†1}, Victoria L Demetriou^{†1}, Maria A Emmanuel^{1,2}, Charalambos K Kakoyiannis² and Leondios G Kostrikis*¹

Address: ¹Department of Biological Sciences, University of Cyprus, Nicosia, Cyprus and ²Department of Veterinary Services, Ministry of Agriculture, Natural Resources and Environment, Nicosia, Cyprus

Email: Andreas V Hadjinicolaou - ah499@cam.ac.uk; Victoria L Demetriou - victoria.demetriou@ucy.ac.cy;

Maria A Emmanuel - emgenic@logos.cy.net; Charalambos K Kakoyiannis - director@vs.moa.gov.cy; Leondios G Kostrikis* - lkostrik@ucy.ac.cy

* Corresponding author †Equal contributors

Published: 19 May 2009

Received: 23 October 2008

BMC Microbiology 2009, 9:97 doi:10.1186/1471-2180-9-97

Accepted: 19 May 2009

This article is available from: <http://www.biomedcentral.com/1471-2180/9/97>

© 2009 Hadjinicolaou et al; licensee BioMed Central Ltd.

This is an Open Access article distributed under the terms of the Creative Commons Attribution License (<http://creativecommons.org/licenses/by/2.0>), which permits unrestricted use, distribution, and reproduction in any medium, provided the original work is properly cited.

Abstract

Background: A fast and simple two-step multiplex real-time PCR assay has been developed to replace the traditional, laborious *Salmonella* serotyping procedure. Molecular beacons were incorporated into the assay as probes for target DNA. Target sequences were regions of the *invA*, *prot6E* and *fliC* genes specific for *Salmonella* spp. *Salmonella* Enteritidis and *Salmonella* Typhimurium, respectively, the two most clinically relevant serotypes. An internal amplification positive control was included in the experiment to ensure the optimal functioning of the PCR and detect possible PCR inhibition. Three sets of primers were used for the amplification of the target sequences. The results were compared to those of the Kauffmann-White antigenic classification scheme.

Results: The assay was 100% sensitive and specific, correctly identifying all 44 *Salmonella* strains, all 21 samples of *S. Enteritidis* and all 17 samples of *S. Typhimurium* tested in this work. Therefore, the entire experiment had specificity and sensitivity of 100%. The detection limit was down to 10 copies of DNA target per 25 µl reaction.

Conclusion: The assay can amplify and analyse a large number of samples in approximately 8 hours, compared to the 4 to 5 days conventional identification takes, and is thus considered a very promising method for detecting the two major serotypes of *Salmonella* quickly and accurately from clinical and environmental samples.

Background

Salmonella is a gram-negative, facultative anaerobic, flagellated bacterium. It is the pathogenic agent of salmonellosis, a major cause of enteric illness and typhoid fever, leading to many hospitalisations and a few rare deaths if no antibiotics are administered. *Salmonella* outbreaks are

linked to unhygienic food preparation, cooking, reheating and storage practices. The bacterium can be isolated from raw meat and poultry products as well as from milk and milk-based products [1]. The detection of *Salmonella* therefore remains a highly important issue in microbiological analysis for food safety and standards.

Because the nomenclature for the *Salmonella* genus is at times confusing, this publication will follow the current literature [2,3]. The CDC [3] distinguishes two *Salmonella* species (or subgenera): *S. enterica* and *S. bongori*. *S. enterica* is further divided into six subspecies, of which *S. enterica* subsp. *enterica* is the most clinically significant, causing 99% of *Salmonella* infections. In the present study we are concerned with its two main serovars: *Salmonella enterica* serovar Typhimurium (group D) denoted *S. Typhimurium*, and *Salmonella enterica* serovar Enteritidis (group B) denoted *S. Enteritidis*, which are the most commonly isolated *Salmonellae* from food-borne outbreaks.

Identification of the disease-causing *Salmonella* serovars is currently a lengthy process, and its initial isolation from food samples can be difficult as the bacteria can be present in small numbers and many closely related bacteria may be found within the same sample [4]. For this reason, pre-enrichment steps are required for all samples [5,6]. The current accepted method for isolation of *Salmonella* from foodstuffs is a well established procedure – ISO 6579, laborious and time-consuming, taking up to 5 days to complete [7,8]. The most widely-used method used to characterise *Salmonella* into its subspecies is the Kauffman-White serotyping system [9], based on the variability of the O, H and Vi antigens [9-11]. Apart from being arduous, this method can not identify a small number of *S. enterica* samples that lack either the O antigen alone or both the O and the H antigens [12]. Therefore there is a need for fast, sensitive and specific "in the field" detection, using nucleic acid-based technologies such as molecular beacon-based real-time PCR, to reduce the time needed to complete the assay, but also improve the level of accuracy and reliability.

In this study, molecular beacons [13-15] and real-time PCR technology are combined to develop a fast, sensitive, clear-cut method of detection of *Salmonella* spp. and its two most clinically significant serotypes from clinical, food and environmental specimens, which reduces the time taken to identify the *Salmonella* strain from an environmental sample, and is precise enough to distinguish between serovars. PCR-based methods targeting various genes are usually more rapid and sensitive than culture-based methods, and the high specificity and high sensitivity of molecular beacons means they can be successfully combined with real-time PCR assays, and so provide a quick, accurate method for detection and analysis, making them ideal for routine diagnosis. Recent studies show that real-time PCR is gradually replacing gel electrophoresis [16-24] as it is suitable for large numbers of samples and involves automatic and fast analysis, as well as being able to execute multiplex protocols. Also, like in all probe-based assays, molecular beacons offer additional specificity. Recent studies have employed molecular beacons in PCR for the detection of *Salmonella* [25-27].

Here the detection of *Salmonella* and the discrimination between *S. Typhimurium* and *S. Enteritidis* serotypes is done by targeting 133–136 nt regions of three genes, while an artificial internal amplification control (IAC) is also incorporated. The target for *Salmonella* spp is *invA*, as it is highly conserved in almost all *Salmonella* serotypes [28,29], and its specificity is apparent from its continuous use in previous studies [18,24,28,30-43]. The target for the specific detection of *S. Enteritidis* is *prot6E*, whose absence from *S. Enteritidis* strains appears to be very rare [18], and the *fliC* gene has been chosen as a single target for the presence of *S. Typhimurium*. The method described here for the detection of *Salmonella* spp. from environmental and food specimens, not only reduces the time taken to identify the *Salmonella* strain, but is also precise enough to distinguish between its clinically significant serovars.

Methods

Bacterial samples

The primary *Salmonella* samples used in this study (Table 1) were obtained from various animal, food and environmental sources at the Cyprus Veterinary Services (Ministry of Agriculture, Natural Resources and Environment, Nicosia, Cyprus), which is the National Reference Laboratory of *Salmonella* for Cyprus. The commercially available bacterial strains listed in Table 2 were obtained from the American Type Culture Collection (ATCC, Manassas, USA), the National Collection of Type Cultures (NCTC, Health Protection Agency, London, UK) and MERCK KGaA (Darmstadt, Germany). The reference *S. enterica* serovars listed in Table 3 were obtained from the Community Reference Laboratory for *Salmonella* at the National Institute for Public Health and the Environment (RIVM, Bilthoven, the Netherlands). Thirty-eight *S. Typhimurium* and *S. Enteritidis* samples as well as six different *Salmonella* serovars have been incorporated to ensure that the assay could correctly identify and differentiate between serotypes of *S. enterica*. Eighteen non-*Salmonella* bacterial samples were also included to ensure accurate experimental design and specificity to the *Salmonella* genus.

Bacterial cultures and serotyping

The detection of *Salmonella* spp. was performed based on the ISO 6579:2002 method. In brief, 25 g of clinical specimen (10 g in the case of minced meat in accordance with the EC regulation 2073/2005 – Microbiological Criteria for Foodstuffs) were added to 225 ml of buffered peptone water in a Stomacher® bag, sealed and placed in a Stomacher® blender for 3 min. The blended sample was incubated for 18 h at 37°C and a 0.1 ml aliquot of sample was inoculated into 10 ml Rappaport-Vassiliadis medium with Soya (RVS) and into 10 ml Muller-Kauffmann tetrathionate/novobiocin (MKTTn) broth; these cultures were incubated for 24 h at 41.5°C and 37°C, respectively. Each culture was inoculated into xylose lysine deoxycho-

Table 1: Primary Bacterial Strains^a

| Bacterial strain | Sample ID | Source of Sample |
|-------------------------------|----------------|-------------------------------------|
| <i>Salmonella</i> Enteritidis | CVS-140/1 | Intestine from beef |
| <i>Salmonella</i> Enteritidis | CVS-141/1-5 | Liver & ovaries from egg layer hens |
| <i>Salmonella</i> Enteritidis | CVS-4054/1 | Lymph ganglions |
| <i>Salmonella</i> Enteritidis | CVS-431/1/1 | Intestine from canaries |
| <i>Salmonella</i> Enteritidis | CVS-4325/4, 5 | Skin from neck of chicken |
| <i>Salmonella</i> Enteritidis | CVS-442/1/1 | Fish food |
| <i>Salmonella</i> Enteritidis | CVS-4516/1 | Veal |
| <i>Salmonella</i> Enteritidis | CVS-4532/1 | Parrot |
| <i>Salmonella</i> Enteritidis | CVS-4540/1 | Parrot |
| <i>Salmonella</i> Enteritidis | CVS-4666/1 | Faeces from egg layer hens |
| <i>Salmonella</i> Enteritidis | CVS-4756/1 | Faeces from hens farmed for meat |
| <i>Salmonella</i> Enteritidis | CVS-4807/1/1-3 | Skin from neck of chicken |
| <i>Salmonella</i> Enteritidis | CVS-4809/2 | Skin from neck of chicken |
| <i>Salmonella</i> Enteritidis | CVS-4980/1 | Faeces from chicken |
| <i>Salmonella</i> Enteritidis | CVS-5212/1 | Faeces from egg layer hens |
| <i>Salmonella</i> Enteritidis | CVS-54/1 | Faeces from egg layer hens |
| <i>Salmonella</i> Enteritidis | CVS-4792/1 | Lymph ganglions |
| <i>Salmonella</i> Enteritidis | CVS-4754/1 | Lymph ganglions |
| <i>Salmonella</i> Enteritidis | CVS-2553/4 | Skin from neck of chicken |
| <i>Salmonella</i> Typhimurium | CVS-3225/1/1-5 | Sheftalia (pork sausage) |
| <i>Salmonella</i> Typhimurium | CVS-4074/1 | Parrot |
| <i>Salmonella</i> Typhimurium | CVS-4076/1 | Pigeon |
| <i>Salmonella</i> Typhimurium | CVS-4255/1 | Beef |
| <i>Salmonella</i> Typhimurium | CVS-4345/4, 5 | Skin from neck of chicken |
| <i>Salmonella</i> Typhimurium | CVS-4979/1 | Dust from egg layer hen cages |
| <i>Salmonella</i> Typhimurium | CVS-4981/1 | Fish meal animal feed |
| <i>Salmonella</i> Typhimurium | CVS-5090/1 | Faeces from finches |
| <i>Salmonella</i> Typhimurium | CVS-55/1 | Faeces from egg layer hens |
| <i>Salmonella</i> Typhimurium | CVS-920/1-3 | Egg yolk |
| <i>Salmonella</i> Typhimurium | CVS-131/2 | Swab from swine |
| <i>Salmonella</i> Typhimurium | CVS-729/2 | Swab from swine |
| <i>Salmonella</i> Typhimurium | CVS-3794/1 | Water |
| <i>Salmonella</i> Typhimurium | CVS-3822/1 | Water |
| <i>Salmonella</i> Typhimurium | CVS-142/1/1 | Lymph ganglions |

^a Identified by culture and serotyping methods as described in the Materials and Methods

late agar (XLD) and brilliant green agar (BGA) and incubated at 37°C for 24 h. One colony was selected from each XLD and BGA plate and spread onto nutrient agar for incubation at 37°C for 24 h. The resulting colonies were subject to biochemical analysis and serotyping. *Salmonella* spp. was characterised into different serovars on the basis of their surface (LPS, O-antigens) and flagellar antigens (H-antigens) as defined by the Kauffman-White Scheme [10,44] and based on the Global Salm-Surv laboratory protocol of the World Health Organisation (Global Salm-Surv, Serotyping of *Salmonella enterica* O and H antigen, Level 3 Training Course, WHO, 6th edition, Jan. 2004).

To extract DNA for use in the molecular detection assay, bacteria were cultured on the XLD agar and one colony was selected and grown on nutrient agar. A colony was then selected and incubated in 5 ml nutrient broth, 1 ml of which was transferred into a 1.5 ml tube for centrifugation for 10 min at 18,000 rcf. The supernatant was dis-

carded and the cell pellet was kept at -80°C until DNA extraction.

Bacterial genomic DNA preparation

Bacterial genomic DNA was extracted from the cell pellets using QIAGEN DNeasy Blood and Tissue Kit (Hilden, Germany) according to the manufacturer's instructions. The purified DNA was eluted in 100 µl of AE buffer and the concentration was determined by measuring the optical density at 260 nm using a NanoDrop UV spectrophotometer (NanoDrop Technologies, USA). The extracted DNA was kept at -30°C until further use.

Internal amplification control

An artificial 129 nt oligonucleotide fragment was designed as an IAC to be amplified by the same primers as the *invA* target. The IAC is a completely synthetic and unique oligonucleotide, designed to avoid sequence homology with any entries in the GenBank database,

Table 2: Commercially Available Strains

| Bacterial Strains | Reference ID |
|-------------------------------|---------------------|
| <i>Salmonella</i> Typhimurium | 14028 ^a |
| <i>Salmonella</i> Enteritidis | 13076 ^a |
| <i>Staphylococcus aureus</i> | 1803 ^b |
| <i>Staphylococcus aureus</i> | 25923 ^a |
| <i>Bacillus cereus</i> | 7464 ^b |
| <i>Bacillus cereus</i> | 11145 ^b |
| <i>Bacillus cereus</i> | 11778 ^a |
| <i>Bacillus subtilis</i> | 110649 ^c |
| <i>Enterobacter aerogenes</i> | 13048 ^a |
| <i>Enterococcus faecalis</i> | 29212 ^a |
| <i>Escherichia coli</i> | 25922 ^a |
| <i>Escherichia coli</i> O157 | 35150 ^a |
| <i>Listeria innocua</i> | 11288 ^b |
| <i>Listeria ivanovie</i> | 11846 ^b |
| <i>Listeria ivanovie</i> | 19119 ^a |
| <i>Listeria monocytogenes</i> | 11994 ^b |
| <i>Micrococcus luteus</i> | 9341 ^a |
| <i>Proteus vulgaris</i> | 13315 ^a |
| <i>Pseudomonas aeruginosa</i> | 27853 ^a |
| <i>Rhodococcus equi</i> | 1621 ^b |

^a Strains obtained from American Type Culture Collection (ATCC), Manassas, USA <http://www.atcc.org>

^b Strains obtained from National Collection of Type Cultures (NCTC), London, UK <http://www.nctc.org.uk>

^c Strains obtained from MERCK KGaA, Darmstadt, Germany <http://www.merck.de>

tested using the BLAST (Basic Local Alignment Search Tool) software [45].

Primer and beacon design

All primers and molecular beacons described here were designed *de novo* for this study (See additional file 1: Oligonucleotide primers and molecular beacons in the real-time PCR assay). Small regions (133–136 nt) of the *invA*, *prot6E* and *fliC* genes were used as target sequences for the detection of *Salmonella* spp. *S. Enteritidis* and *S. Typhimurium*, respectively. The primers and the molecular beacons were designed based on sequences of the above genes found in the GenBank database [http://](http://www.ncbi.nlm.nih.gov/Genbank/index.html)

Table 3: *Salmonella enterica* serovars^a

| Serovar | Reference ID |
|-------------------------------|--------------|
| <i>Salmonella</i> Bredeney | 1030/1 |
| <i>Salmonella</i> Infantis | 1030/4 |
| <i>Salmonella</i> Anatum | 1030/5 |
| <i>Salmonella</i> Hadar | 1030/6 |
| <i>Salmonella</i> Newport | 1030/7 |
| <i>Salmonella</i> Typhimurium | 1030/10 |
| <i>Salmonella</i> Virchow | 1030/11 |
| <i>Salmonella</i> Enteritidis | 1030/17 |

^a Obtained from the Community Reference Laboratory for *Salmonella*, RIVM, Bilthoven, the Netherlands <http://www.rivm.nl>

www.ncbi.nlm.nih.gov/Genbank/index.html using BLAST [45]. The molecular beacons, target oligonucleotides and primers were synthesised by MWG-Biotech AG Ltd (Ebersberg, Germany) and the Midland Certified Reagent Company, Incorporated (Texas, USA). For the stem formation, the ends of each beacon were designed to have a high GC content and to be complementary to each other. All beacons were labelled with DABCYL, i.e., 4'-(4'-dimethylaminophenylazo) benzoic acid at the 3' end and with one of four fluorophores at the 5' end. Molecular beacon MB_{invA}, was labelled with the fluorophore FAM (Fluorescein); MB_{prot6E} with TET (Tetrachloro-6'-carboxyfluorescein); MB_{fliC} with HEX (hexachlorofluorescein); and MB_{IAC}, with ROX (6'-carboxy-X-rhodamine). Within the loop element, each beacon contains a 22–25 nucleotide-long probe sequence complementary to the target. In addition to the probe sequence, each beacon has 4–5 of the 12 bases of its two arms also complementary to the target. In this non-traditional way, once the beacon is in the open structure, it binds more forcefully to the target and the hybrid formed is more stable, and the maximum distance possible between fluorophore and quencher is created.

Thermal denaturation characteristics of the molecular beacons

To assess the thermodynamic characteristics, the quality and the purity of the molecular beacons used in this study, a melting curve analysis was performed on each using the 7900 HT Fast Real-Time PCR System (Applied Biosystems, Foster City, CA, USA). Briefly, the reaction consisted of a 25 µl solution containing 12.5 µl Platinum[®] PCR Supermix (Invitrogen), 1 µl of the beacon probe at the appropriate concentration with or without 1 µl (100 pmol) of a single-stranded oligonucleotide perfectly complementary to the probe. The cycling parameters were as follows: 1 cycle for 2 min at 95°C followed by 50 cycles, each consisting of the data collection step for 30 s and a second step for 10 s, each starting at 80°C and employing auto-incrementation of -1°C per half-minute cycle until 31°C. Changes in fluorescence were measured at 490 nm and the data was collected at each temperature interval.

PCR target standards synthesis, amplification and quantification

PCR target standards of the *fliC*, *invA*, *prot6E* and IAC target sequences were synthesised by PCR amplification using long overlapping primers. Each PCR reaction was performed in a 50-µl reaction volume containing 45 µl of Platinum[®] PCR Supermix (Invitrogen, Inc., Carlsbad, CA, USA) and 25 pmoles each of the following primer pairs: for the *fliC* target, TF_{fliC} and TR_{fliC}; for the *invA* target, TF_{invA} and TR_{invA}; for the *prot6E* target, TF_{prot6E} and TR_{prot6E}; and for the IAC, TFIAC and TRIAC (See additional file 1: Oligonucleotide primers and molecular bea-

cons in the real-time PCR assay). Amplification was performed with an activation step of 94°C for 30 s, followed by 20 cycles, each consisting of 94°C for 20 s, 68°C for 30 s and 72°C for 20 s, followed by a final extension step of 72°C for 5 min in an Eppendorf Mastercycler (Eppendorf AG, Hamburg, Germany). Three µl of the product from the first PCR was used in a secondary PCR in a 50-µl reaction volume containing 1 × of Platinum® PCR Supermix (Invitrogen, Inc., Carlsbad, CA) and 20 pmoles each of the following PCR primer pairs: for the *fliC* target, 585 and 717; for the *invA* target, 302 and 437; for the *prot6E* target, 438 and 572; and for the IAC, 302 and 437 (See additional file 1: Oligonucleotide primers and molecular beacons in the real-time PCR assay). Amplification was performed with an activation step of 94°C for 30 s, followed by 40 cycles, each consisting of 94°C for 20 s, the annealing temperature for 30 s and 72°C for 20 s, followed by a final extension step of 72°C for 5 min in an Eppendorf Mastercycler (Eppendorf AG, Hamburg, Germany). The annealing temperature for the *fliC* primers was 59°C, for the *invA* primers 58°C, for the *prot6E* primers 56°C and for the IAC primers 58°C. The resulting product was then cleaned using the QIAquick PCR Purification Kit (QIAGEN GmbH, Hilden, Germany) and eluted in 50 µl of EB buffer. The PCR products were run on a 2% agarose gel with a 50 bp DNA ladder (Invitrogen) and the DNA concentration of each was measured on the NanoDrop ND-1000 UV Spectrophotometer (Wilmington, DE). The number of molecules per unit volume was calculated from the measured concentration and the molecular weight of each oligonucleotide. The amplified targets were then diluted to concentrations of 10⁶, 10⁵, 10⁴, 10³, 10² and 10 copies per 5 µl to be used as target standards of known concentration.

Standard curves

Uniplex real-time PCR reactions were performed on 10-fold serial dilutions of the PCR targets, synthesised and prepared as described above. Reactions of 25 µl volume were set up, containing 12.5 µl Platinum® qPCR Supermix-UDG (Invitrogen, Carlsbad, CA), 1 µl of forward primer and 1 µl of reverse primer (20 pmol/µl), 1 µl of the corresponding molecular beacon at the concentration determined appropriate from the melting curve analysis (4.9 pmol/µl MBinVA, 10 pmol/µl MBfliC, 4.4 pmol/µl MBprot6E and 50 pmol/µl MBIAC) 4.5 µl H₂O and 5 µl of the PCR target standard. The reactions were performed on the 7900 HT Real-Time PCR System (Applied Biosystems, Foster City, CA, USA) and the cycling parameters were as follows: activation step at 95°C for 10 min, followed by 50 cycles, each consisting of 95°C for 15 s, 50°C for 30 s (data collection step) and 72°C for 30 s.

Uniplex real-time PCR

The real-time PCR analysis was made with by the 7900 HT Fast Real-Time PCR System (Applied Biosystems) using

the Platinum® Quantitative PCR SuperMix-UDG (Invitrogen) on all of the samples described above. Each 25 µl uniplex PCR reaction contained 5 µl of the extracted DNA, and was carried out as described above. The fluorescence given out on hybridisation between each beacon and its target DNA was measured directly and the resulting amplification curves were processed immediately with the 7900 HT Sequence Detection Systems software v2.2.2 (Applied Biosystems, Foster City, CA). To verify that the fluorescence signals were due to PCR amplification of the template DNA and not any other contaminant, negative or non-template controls were also run, where sterile water replaced the DNA template in the reaction mixture.

Double duplex real-time PCR

Having tested all sets of beacons and primers in uniplex reactions, the samples were run again in a two-step duplex assay. In step 1, 25 µl reactions were set up, containing 12.5 µl of Platinum Quantitative Supermix-UDG (Invitrogen), 1 µl of each of primers 302 and 437 (20 pmol/µl), 1 µl of MBIAC (50 pmol/µl), 1 µl of MBinVA (4.9 pmol/µl), 0.5 µl of the synthetic IAC (2 × 10⁵ copies/µl). To this, 2 µl of 100-fold dilution of sample DNA were added and the volume was made up with sterile water or, in the case of non-template controls, the sample DNA was replaced with sterile water. In step 2, each reaction had a total volume of 25 µl consisting of 12.5 µl of Platinum Quantitative Supermix-UDG (Invitrogen), 1 µl of each of 572, 585 and 717 (20 pmol/µl), 1 µl of MBprot6E (4.4 pmol/µl) and 2 µl of MBfliC (10 pmol/µl). The final volume was reached by the addition of 2 µl of sample DNA and 3.5 µl of sterile water or, in the case of non-template negative control reactions, 5.5 µl of sterile water only. For both steps, PCR cycling conditions were as described for the standard curve analysis and uniplex reactions. The fluorescence given out on hybridisation between beacon and its target was measured at each cycle.

Results

Thermal denaturation characteristics of molecular beacons

Normalised fluorescence signals for both the beacon and the beacon-target hybrid were plotted against temperature to give a thermal denaturation profile for each beacon (Fig. 1). These profiles were created using an ABI 7900 HT Fast Real-Time PCR System (Applied Biosystems, Foster City, CA) to determine the optimal hybridisation temperature between the beacon and its target sequence. Perfectly complementary beacon-target hybrids exist at lower temperatures giving out a bright fluorescence signal. A progressive increase in temperature causes the hybrids to dissociate, followed by a marked decrease in fluorescence. Conversely, the beacons alone unravelled at high temperatures and exhibited a melting temperature above 60°C in all cases. In the temperature interval from 31 to 55°C, the probe-target hybrids elicited significantly stronger fluores-

cence than the probe alone, thus allowing the detection of target sequence at these temperatures. A temperature of 50°C was chosen as an optimal annealing temperature for subsequent real-time PCR studies. At this temperature the difference in fluorescence signal between beacon alone and beacon-target hybrids is large; in the absence of target any fluorescence detected is background level and the temperature is high enough to prevent less energetically favourable hybrids from forming, e.g., primer dimers or

beacon-primer dimers. In the process of carrying out the melting curve analysis for all beacons, different concentrations were tested, to find the appropriate concentration at which the fluorescence signal was neither too low nor saturated. The concentrations at which the particular beacons exhibited the desired amount of fluorescence signal in these reactions were: MBIAC, 50 pmol/μl; MBinvA, 4.9 pmol/μl; MBprot6E, 4.4 pmol/μl; and MBfliC, 10 pmol/μl. Finally, these thermal denaturation profiles illustrate

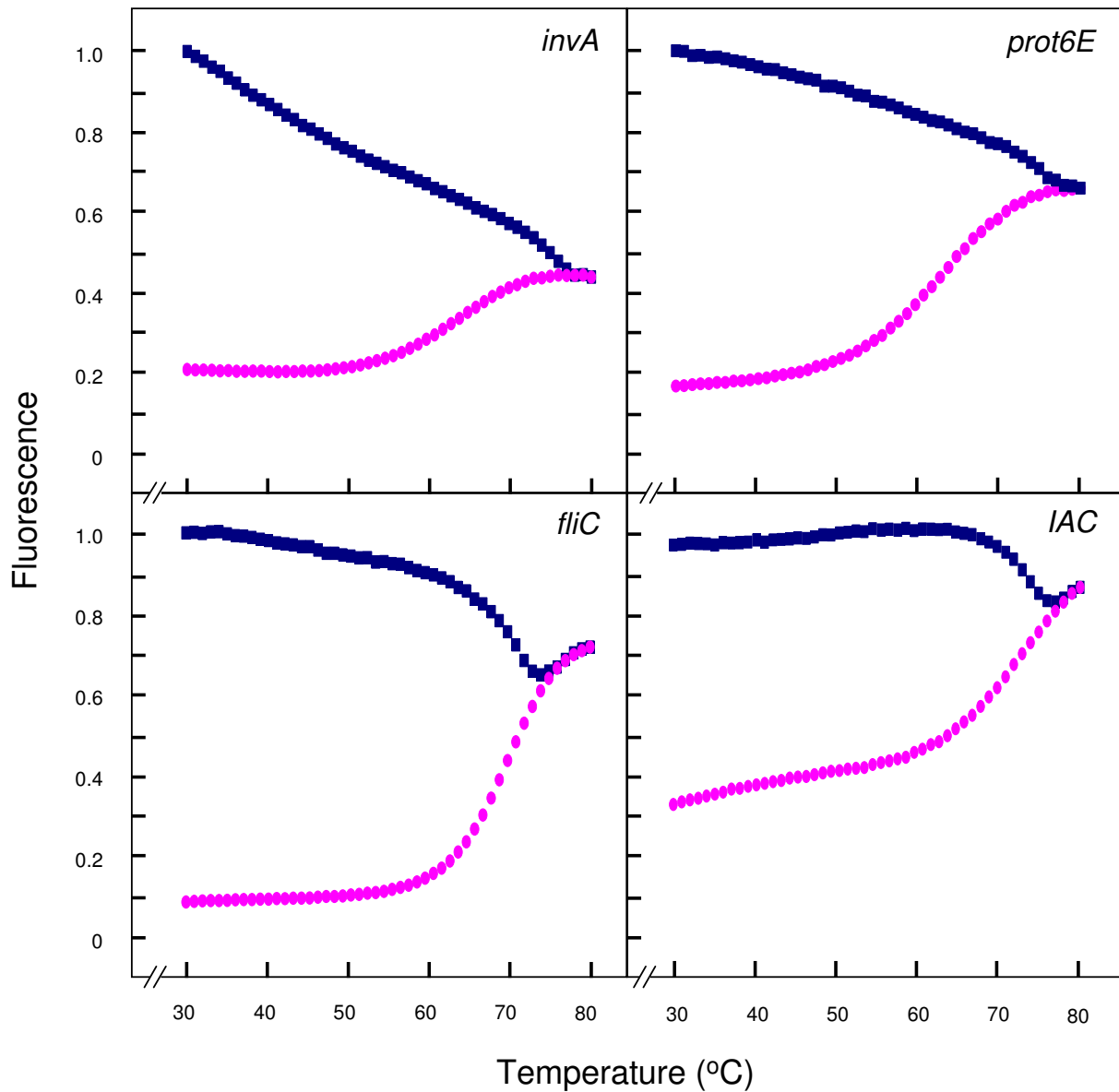


Figure 1
Thermal denaturation profiles of the molecular beacons. Thermal denaturation profiles of the molecular beacons used in this study as established by melting curve analysis (described in Materials and Methods). The figure shows normalised fluorescence thermal transitions of molecular beacon plotted in pink circles and beacon-target complexes plotted in blue squares.

the good quality of the molecular beacons and their efficiency in hybridising with the appropriate target sequence.

Standard curves and limit of detection

Standard curves were initially plotted to ensure the ability of each molecular beacon to detect its specific *Salmonella* target and the detection limit of the assay. The copy numbers of target standards used ranged from 10^1 to 10^6 copies per reaction. These plots represent how the amplification of DNA progresses with each log increase of target copy number. The small standard errors calculated from multiple values of the threshold cycle at which significant DNA amplification was observed (threshold cycle, C_T) for each reaction and indicated on the graphs with horizontal lines above and below each plotted point, suggest that the PCR amplification is highly reproducible. The C_T values for the target sequences depended on the initial DNA amount in each reaction as shown by the linear relationship of standard curves along a 6-log range which yields an R^2 correlation value higher than 0.994 in all three cases (Fig. 2). The correlation was 0.995 with 76% efficiency for *invA*, 0.997 and 84% efficiency for *prot6E* and 0.999 and 100% efficiency for *fliC*. As the reactions worked well for all target standard concentrations tested, the lower limit of detection for the assay was set to be 10 copies of the required target fragment per reaction. Based on the standard curves and the limit of detection of this assay, negative results were defined as those exhibiting C_T values higher than 45.

Detection of *S. enterica* alleles in bacterial samples by molecular beacon-based uniplex real-time PCR

The molecular beacon-based real-time PCR assay designed in this study was tested on environmental and food samples of *S. Enteritidis* and *S. Typhimurium* (Table 1), as well as several commercially available bacterial strains (Table 2) and various *Salmonella* serovars obtained from a reference laboratory for *Salmonella* (Table 3). All samples were investigated first by uniplex assays to detect *invA*, *prot6E* and *fliC* (Table 4). In the reaction for detection of *invA*, all 44 *Salmonella* samples were positive and all 18 non-*Salmonella* samples were undetectable. Positive results (≤ 10 copies of DNA per reaction) had C_T values ranging from 15 to 25. In the *prot6E* reaction, all 21 *S. Enteritidis* samples gave positive PCR results and all 41 non-*Enteritidis* samples were negative. Positive samples for the *prot6E* gene had C_T values ranging between 15 and 18 with one exception, the commercially available specimen of *S. Enteritidis* (Table 3) for which fluorescence detection significantly increased around cycle 30. Finally, in the *fliC* reaction, all 17 *S. Typhimurium* samples gave positive PCR results and all 45 non-*Typhimurium* samples were negative. Positive results had C_T values ranging from 15 to 18 cycles. These results showed that the primers and beacons for each reaction work well individually and that they amplify and detect their target sequence with very high specificity and sensitivity. The C_T values exhibited by the samples in these experiments, compared to the plot of the standards of known concentration, indicated that the extracted DNA from the bacterial samples

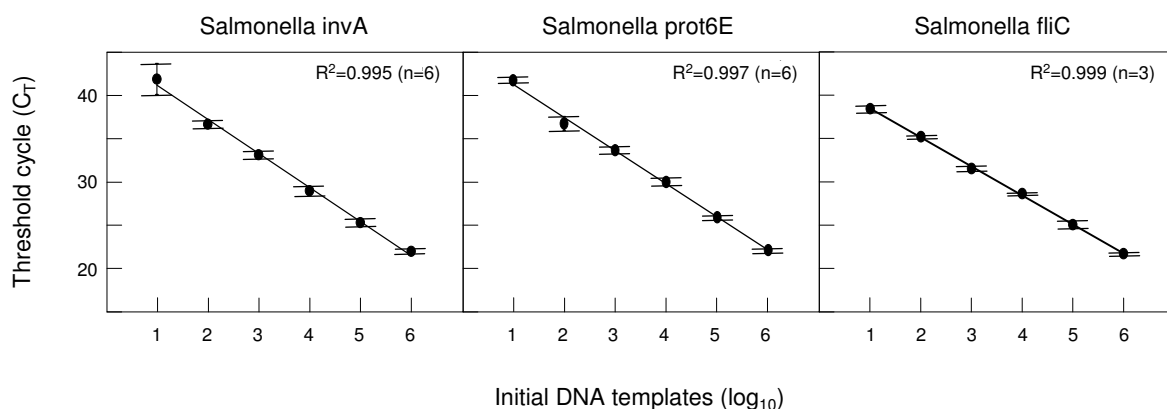


Figure 2

Standard curves for targets *invA*, *fliC* and *prot6E*. Standard curves for molecular beacon-based real-time PCR detection of targets *invA*, *fliC* and *prot6E*. The plots illustrate the relationship of known number of target DNA copies per reaction to the threshold cycle of detection (C_T) for each molecular beacon reaction. The C_T is directly proportional to the log of the input copy equivalents, as demonstrated by the standard curves generated.

Table 4: Real-time PCR results obtained using the available bacterial strains

| Bacterial strain ^a | Reference ID | PCR Result ^b | | |
|-------------------------------|---------------|-------------------------|-------------|---------------------------|
| | | IAC ^c | <i>invA</i> | <i>prot6E</i> <i>fliC</i> |
| <i>Salmonella</i> Enteritidis | CVS-140/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-141/1-5 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4054/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4311/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4325/4, 5 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4421/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4516/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4532/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4540/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4666/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4756/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4807/1-3 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4809/2 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4980/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-5212/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-54/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4792/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-4754/1 | + | + | - |
| <i>Salmonella</i> Enteritidis | CVS-2553/4 | + | + | - |
| <i>Salmonella</i> Typhimurium | CVS-3225/1-5 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-4074/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-4076/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-4255/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-4345/4, 5 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-4979/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-4981/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-5090/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-55/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-920/1-3 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-131/2 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-729/2 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-3794/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-3822/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | CVS-1421/1 | + | - | + |
| <i>Salmonella</i> Typhimurium | 14028 | + | - | + |
| <i>Salmonella</i> Enteritidis | 13076 | + | + | - |
| <i>Staphylococcus aureus</i> | 1803 | - | - | - |
| <i>Staphylococcus aureus</i> | 25923 | - | - | - |
| <i>Bacillus cereus</i> | 7464 | - | - | - |
| <i>Bacillus cereus</i> | 11145 | - | - | - |
| <i>Bacillus cereus</i> | 11778 | - | - | - |
| <i>Bacillus subtilis</i> | 110649 | - | - | - |
| <i>Enterobacter aerogenes</i> | 13048 | - | - | - |
| <i>Enterococcus faecalis</i> | 29212 | - | - | - |
| <i>Escherichia coli</i> | 25922 | - | - | - |
| <i>Escherichia coli</i> O157 | 35150 | - | - | - |
| <i>Listeria innocua</i> | 11288 | - | - | - |
| <i>Listeria ivanovie</i> | 11846 | - | - | - |
| <i>Listeria ivanovie</i> | 19119 | - | - | - |
| <i>Listeria monocytogenes</i> | 11994 | - | - | - |
| <i>Micrococcus luteus</i> | 9341 | - | - | - |
| <i>Proteus vulgaris</i> | 13315 | - | - | - |
| <i>Pseudomonas aeruginosa</i> | 27853 | - | - | - |
| <i>Rhodococcus equi</i> | 1621 | - | - | - |
| <i>Salmonella</i> Bredeney | 1030/1 | + | - | - |
| <i>Salmonella</i> Infantis | 1030/4 | + | - | - |

was higher than the range of concentrations tested by the standards ($>10^7$ copies per reaction). Therefore 100-fold dilutions of all extracted DNA samples were prepared for use in the two-step duplex assay, so that the resulting C_T values would fall within the range seen on the standard curves.

Detection of *S. enterica* alleles in bacterial samples by molecular beacon-based two-step duplex real-time PCR

All samples were subsequently tested employing a double duplex strategy, where the first step of the assay involved a real-time PCR reaction containing the primers and beacons for the amplification and detection of both the *invA* gene and the IAC, a quantity of which was inserted into each reaction. This step of the assay determined the presence or absence of the *invA* gene to verify whether the sample was *Salmonella* or not. Simultaneous amplification and detection of the IAC in all samples ensured accurate PCR performance, therefore excluding the possibility of false-negative results. Amplification and detection of both *invA* and the IAC were clear in all *Salmonella* samples, whereas only the IAC amplification was detected in non-*Salmonella* samples. Representative amplification plots from *Salmonella* and other bacteria for the first step reaction are seen in Fig. 3. The results demonstrate that this reaction correctly recognises samples in which *Salmonella* exist from samples in which it does not.

All samples found positive for *invA* in the first step were then tested in the second step of the assay, another duplex real-time PCR reaction containing the components for amplification and detection of both *prot6E* and *fliC* targets. In all *S. Typhimurium* samples *fliC* was the only target detected, in all *S. Enteritidis* samples *prot6E* was the only target detected and in all other *Salmonella* samples, both targets were undetected. The results show that this reaction clearly and accurately distinguishes between *S. Typhimurium* strains, *S. Enteritidis* strains and other *Salmonella* serotypes. Representative amplification plots from *S. Typhimurium*, *S. Enteritidis* and other *Salmonellae* for the second step reaction are seen in Fig. 4, clearly showing that the *prot6E* and *fliC* components designed in this study work well together in a multiplex real-time PCR reaction.

In both the uniplex and double duplex assays, non-template controls were included to verify the absence of false-positive results. In all cases they exhibited undetectable amplification of the targets ($C_T > 45$).

Selectivity of the real-time assay

The selectivity and accuracy of the test is measured by calculating the values for specificity and sensitivity. Specificity is the probability that the PCR will be negative among specimens that should not possess the gene and is calcu-

Table 4: Real-time PCR results obtained using the available bacterial strains (Continued)

| | | | | |
|-------------------------------|---------|---|---|---|
| <i>Salmonella</i> Anatum | 1030/5 | + | - | - |
| <i>Salmonella</i> Hadar | 1030/6 | + | - | - |
| <i>Salmonella</i> Newport | 1030/7 | + | - | - |
| <i>Salmonella</i> Typhimurium | 1030/10 | + | - | + |
| <i>Salmonella</i> Virchow | 1030/11 | + | - | - |
| <i>Salmonella</i> Enteritidis | 1030/17 | + | + | - |

^a The origin of these samples has been described in tables 1, 2 and 3.

^b Detectable real-time PCR amplification signal is denoted by the symbol (+) which indicates C_T values ≤ 45 . Undetectable signal (C_T values >45) is denoted by (-) and indicates the absence of template DNA.

^c IAC is an artificial internal amplification control that does not correspond to any of the submitted sequences in the Genbank.

lated using the formula: true negative/(true negative + false positive). Sensitivity shows the strength of the test in recognising what we are looking for, i.e., in correctly identifying the specific serotype. The formula used for estimation of sensitivity is: true positive/(true positive + false negative). For the reaction targeting the *invA* gene, all 44 *Salmonella* samples investigated were positive indicating a sensitivity of 100%. The specificity was also 100% since all non-*Salmonella* samples gave negative results, with undetectable fluorescence signals after 50 cycles. In the *prot6E* reaction all *S. Enteritidis* samples analysed were identified correctly with positive PCR results and all non-*Enteritidis* samples were negative for this target. Thus, this reaction also had sensitivity and specificity of 100%. Similarly, in the assay for *fliC* detection, all *S. Typhimurium* samples tested were positive for the target. The assay's sensitivity was 100%, matched by an equal specificity value as all non-*Typhimurium* samples tested gave negative PCR results.

Discussion

Traditional serotyping of *S. enterica* is based on the detection of certain antigens using microbiological techniques and culturing, which are time-consuming and laborious. This study exploits real-time PCR, molecular beacons and genetic variation between different serotypes to devise a quick, accurate and simple assay to reliably identify a bacterial sample as *Salmonella enterica* and further distinguish it as serotypes *S. Typhimurium* or *S. Enteritidis*, the two serovars most commonly associated with food-borne gastroenteritis. The assay described in this study can analyse a large number of samples very quickly, and can also identify as few as 10 copies of target DNA per reaction, potentially even in the presence of thousands of copies of other serotypes. The short completion time of this assay and the ability of PCR to detect even dead bacterial cells, even though these are no longer infectious, highlight the drawbacks of the culturing method and the need to start using modern molecular techniques as opposed to conventional microbiology methods where applicable. PCR sensitivity is superior to that of the bacteriological culturing

methods, as it can detect *Salmonellas* with atypical biochemical features, reducing false-negative results, and it will not mistakenly detect non-*Salmonella* bacteria, reducing the chances of false-positive data [27]. However, further research is necessary to ensure that molecular assays alone can efficiently detect *Salmonella* spp. and its serotypes.

A variety of bacterial samples were used to test the specificity of the assay in the detection of the genus *Salmonella*. At the same time a number of *Salmonella* strains were included to ensure that the detection tests for serovars *S. Typhimurium* and *S. Enteritidis* were specific. The study includes strains from clinical and environmental samples as well as commercially available strains, and a significant number of *S. Typhimurium* and *S. Enteritidis* samples as well as other *Salmonella* serotypes and non-*Salmonella* bacteria. This broad range of samples was included to test the efficacy of the assay. Three genes were specifically targeted: the *invA* gene specific to the genus *Salmonella*; the *prot6E* gene specific to *S. Enteritidis*; and the *fliC* gene specific to *S. Typhimurium*. Due to its specificity, the *invA* gene is an excellent potential target for the detection of *S. enterica* strains alone [18,24,28,30-43]. The fact that it is unique for *S. enterica* and rarely absent from it [46], ensures high specificity and sensitivity in detection assays. The *prot6E* gene is located on a highly conserved, low copy number, 60-kb virulence plasmid specific to *S. Enteritidis* and its absence appears to be very rare [18]. Finally, the *fliC* gene codes for the H1 antigen of *Salmonella*. Targeting the *fliC-i* allele greatly increases the specificity for *S. Typhimurium* identification. In order to detect *S. Typhimurium* with the highest specificity, three genes could ideally be targeted, coding for antigens O:4, H1:i and H2:1,2, as it is the only serovar with this antigen combination [47]. However, this would not only raise the costs of the assay but would compromise the simplicity of design and the potential for including further molecular beacons in the multiplex reaction for identification of other target serotypes. Thus, in this study, as in some other studies [48,49], the *fliC* gene has been chosen as a single target for the presence of *S. Typhimurium*. By designing the *fliC* beacon using a *S. Typhimurium* sequence from the GenBank database as a template, the assay exhibits high sensitivity. However, although it performed with 100% specificity with this particular set of samples, in a different set of samples, e.g., with other *S. enterica* serotypes like Kentucky that also possess the H1:i antigen, it might yield different results. Many different genes have been targeted in previous studies [16,22,25,26,30,47,48,50-54]. However, the above targets did not prove to be specific enough for unique detection and identification.

The IAC used is a synthetic and unique oligonucleotide designed *de novo* for this study. The fact that this IAC is co-

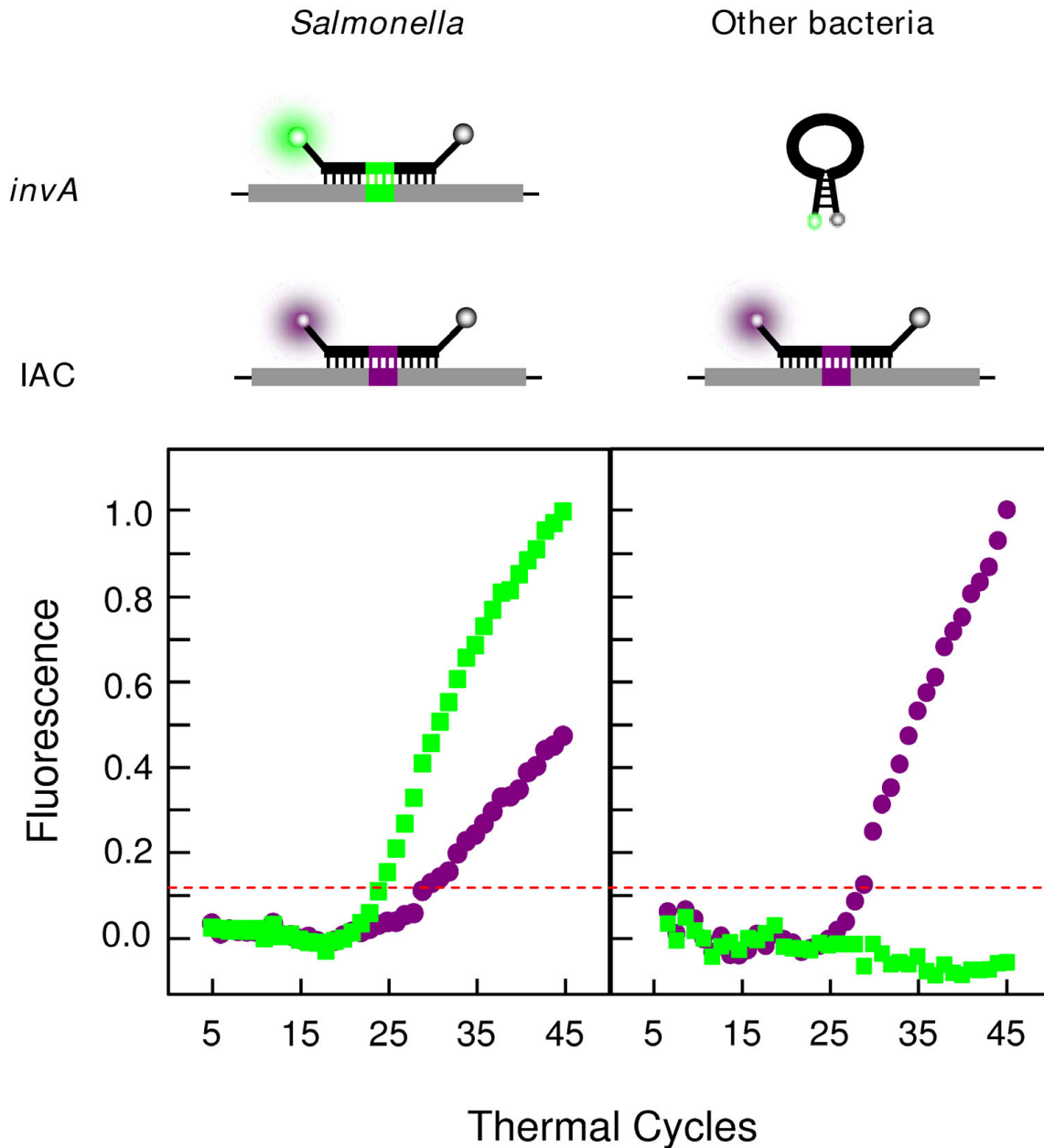
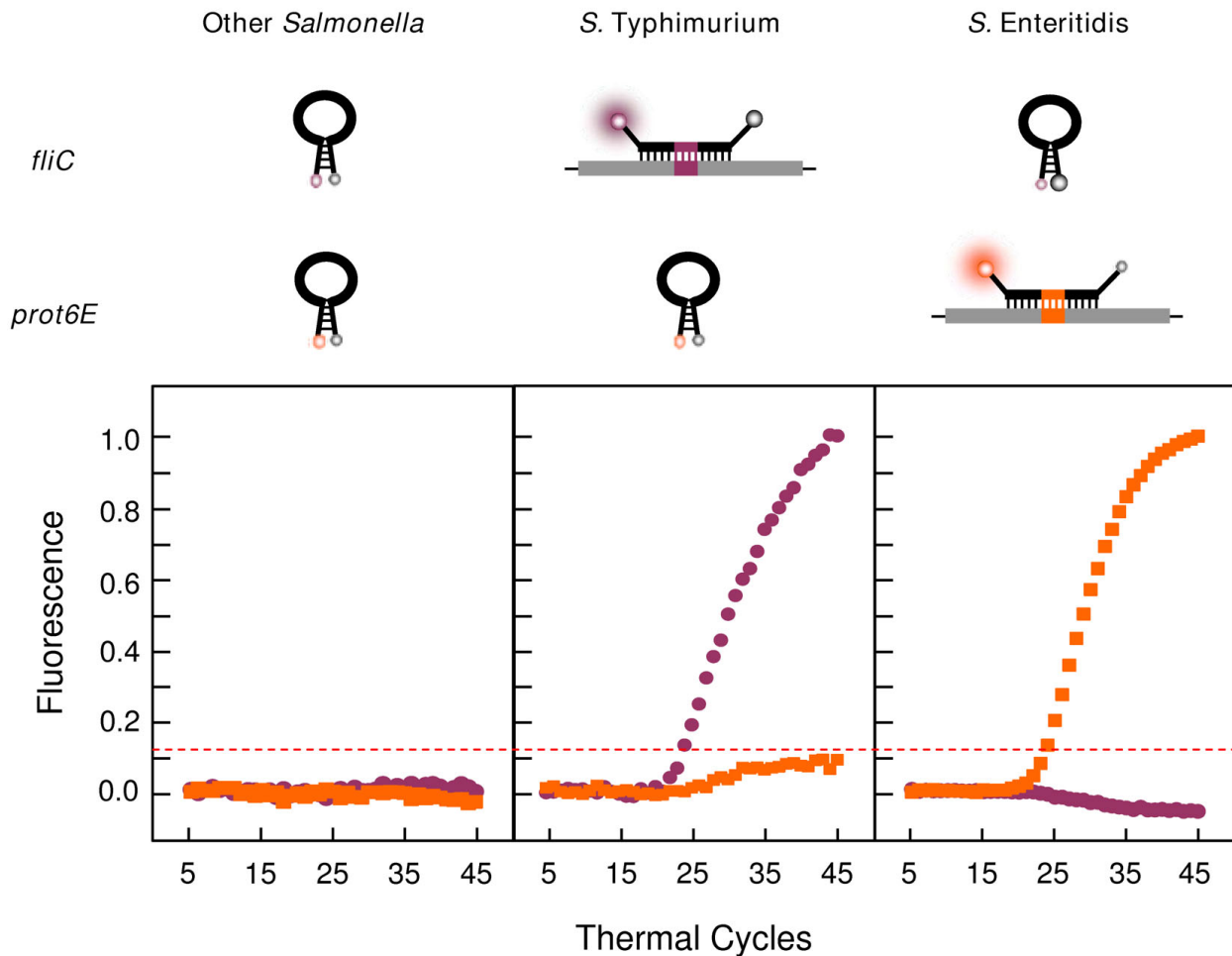


Figure 3
Schematic real-time PCR results for the first step reaction. Representative real-time PCR results as established by the first step multiplex reaction (described in Materials and Methods). The plots show average normalised linear amplification of representative samples shown for demonstration of typical results obtained from *Salmonella* and non-*Salmonella* bacteria. With DNA from non-*Salmonella* bacterial samples, only the IAC-specific, ROX-labelled molecular beacons hybridise to the IAC amplicons, generating violet fluorescence, whereas the *invA*-specific, FAM-labelled molecular beacons retain their stem-and-loop structure and cannot produce a green fluorescent signal. With DNA from *Salmonella* samples, both molecular beacons hybridise to their respective target amplicons and generate both green and violet fluorescence. The dashed line on the plots represents the normalised threshold for detection of fluorescence, the baseline above which fluorescence increases significantly on amplification and detection of the target sequence.

**Figure 4**

Schematic real-time PCR results for the second step reaction. Representative real-time PCR results as established by the second step multiplex reaction (described in Materials and Methods). The plots show average normalised linear amplification of representative samples shown for demonstration of typical results obtained from *S. Typhimurium*, *S. Enteritidis* and other *Salmonella* samples. With DNA from *S. Typhimurium* strains, only *fliC*-specific, HEX-labelled molecular beacons hybridise to the amplicons, generating pink fluorescence, whereas the *prot6E*-specific, TET-labelled molecular beacons retain their stem-and-loop structure and cannot produce an orange fluorescent signal. With DNA from *S. Enteritidis* strains, the *prot6E*-specific, TET-labelled molecular beacons hybridise to their target amplicons and produce an orange fluorescent signal, whereas the *fliC*-specific, HEX-labelled molecular beacons remain dark. With DNA from other *Salmonella* serotypes, no target amplicons are detected and both molecular beacons remain dark. The dashed line on the plots represents the normalised threshold for detection of fluorescence, the baseline above which fluorescence increases significantly on amplification and detection of the target sequence.

amplified with the *invA* fragment using the same primer set but detected by a distinct beacon, does not appear to alter the precision and accuracy of the real-time PCR, and quantification of original target DNA is still possible even in the presence of the control. However, the standard curve protocol for *invA* in the presence of the IAC should be performed for a correct quantitative approach if the assay is to be used for quantification. The *invA* gene has

been used as an internal amplification control in other studies [18], but its application is limited to *Salmonella* assays alone. Furthermore, it has been found that in some rare cases, this gene may be absent and is therefore unreliable as an amplification control even for studies incorporating *Salmonella* specimens alone. This IAC sequence matches no organism in the NCBI libraries and could potentially be used in any such detection assays.

The assays for the *invA*, *fliC* and *prot6E* genes all had a sensitivity and specificity score of 100%. All 45 *Salmonella* samples were positive, with 100% sensitivity. Positive results (>10 copies of DNA per reaction) had C_T values ranging from 15 to 25. One exception, the commercially available specimen of *S. Enteritidis* (Table 3), had a C_T value of approximately 30. Since the *prot6E* gene is located on a virulence plasmid, its absence would not be surprising. Plasmid profiling should be performed to explain the unusually high C_T value observed for this specific specimen. This raises the question of whether selecting a target on a plasmid is a wise choice, but this absence of this plasmid has been found to be rare from *S. Enteritidis* and the low copy numbers (1–2) of the plasmid in the cells make possible the conversion of the assay to a quantitative one which would be correct to a factor of 2. Therefore, using this target for quantification would depend on the accuracy required.

Our study is the first to incorporate four molecular beacons with real-time PCR in a double duplex PCR protocol to detect *Salmonella* spp., *Salmonella* Enteritidis and *Salmonella* Typhimurium in a single assay. Strong fluorescence signals were observed in all positive PCR results in both the uniplex and the duplex assays, indicating the efficiency of the design in the primers and beacons. The sensitivity and specificity of the design and procedure described here give the assay the potential to be converted into a quantitative method, directly applied to samples without the requirement of pre-enrichment stages, making use of the standard curves.

Conclusion

This molecular beacon-based real-time PCR assay is extremely useful for any laboratory in possession of a real-time PCR. It is a fast, reproducible, simple, specific and sensitive way to detect nucleic acids, which could be used in clinical diagnostic tests in the future. The design of the assay gives it potential to be used for quantification, for detection of multiple other serotypes of *Salmonella* or to be modified for the detection of other bacterial samples. Also, more significantly, the sensitivity of the test and its confirmed low limit of detection, are promising factors for the important switch to direct detection from real clinical and environmental samples which have not been previously cultured and have low numbers of bacteria.

Authors' contributions

AVH participated in the assay design, sample preparation, real-time PCR experimental procedures, the analysis and interpretation of the results and drafted the manuscript. VLD carried out sample preparation, real-time experimental procedures, analysis and interpretation of results and drafted the manuscript. MAE carried out the bacterial culturing and serotyping techniques, sample selection, bacte-

rial pellets isolation and helped with the manuscript preparation. CKK participated in sample selection and donated samples for this study. LGK conceived and designed the assay, coordinated the study and participated in sample selection and analysis and interpretation of results. All authors read and approved the final manuscript.

Additional material

Additional File 1

Oligonucleotide primers and molecular beacons in the real-time PCR assay. Table of primer and molecular beacon sequences used in this study.

Click here for file

[<http://www.biomedcentral.com/content/supplementary/1471-2180-9-97-S1.doc>]

Acknowledgements

The authors would like to thank the staff at the Department of Veterinary Services, Nicosia for assisting in sample collection; S. Gilliland and J. Hezka for technical assistance.

This work was supported by research funds from the University of Cyprus (awarded to L. G. Kostrikis) and the Birch Biomedical Research LLC.

References

- Gorman R, Adley CC: **Characterization of *Salmonella enterica* serotype Typhimurium isolates from human, food, and animal sources in the Republic of Ireland.** *J Clin Microbiol* 2004, **42(5)**:2314-2316.
- Tindall BJ, Grimont PA, Garrity GM, Euzéby JP: **Nomenclature and taxonomy of the genus *Salmonella*.** *Int J Syst Evol Microbiol* 2005, **55(Pt 1)**:521-524.
- Brenner FW, Villar RG, Angulo FJ, Tauxe R, Swaminathan B: ***Salmonella* nomenclature.** *J Clin Microbiol* 2000, **38(7)**:2465-2467.
- Baylis CL, MacPhee S, Betts RP: **Comparison of methods for the recovery and detection of low levels of injured *Salmonella* in ice cream and milk powder.** *Lett Appl Microbiol* 2000, **30(4)**:320-324.
- Baylis CL, MacPhee S, Betts RP: **Comparison of two commercial preparations of buffered peptone water for the recovery and growth of *Salmonella* bacteria from foods.** *J Appl Microbiol* 2000, **89(3)**:501-510.
- Hoorfar J, Baggesen DL: **Importance of pre-enrichment media for isolation of *Salmonella* spp. from swine and poultry.** *FEMS Microbiol Lett* 1998, **169(1)**:125-130.
- Uyttendaele M, Vanwildemeersch K, Debevere J: **Evaluation of real-time PCR vs automated ELISA and a conventional culture method using a semi-solid medium for detection of *Salmonella*.** *Lett Appl Microbiol* 2003, **37(5)**:386-391.
- Voogt N, Raes M, Wannet WJ, Henken AM, Giessen AW van de: **Comparison of selective enrichment media for the detection of *Salmonella* in poultry faeces.** *Lett Appl Microbiol* 2001, **32(2)**:89-92.
- Popoff MY, Le Minor L: **Antigenic formulas of the *Salmonella* serovars, 7th revision.** Institut Pasteur, Paris; 1997.
- Popoff MY, Bockemuhl J, Gheesling LL: **Supplement 2002 (no. 46) to the Kauffmann-White scheme.** *Res Microbiol* 2004, **155(7)**:568-570.
- Yoshida C, Franklin K, Konczyk P, McQuiston JR, Fields PI, Nash JH, Taboada EN, Rahn K: **Methodologies towards the development of an oligonucleotide microarray for determination of *Salmonella* serotypes.** *J Microbiol Methods* 2007, **70(2)**:261-271.

12. Hoorfar J, Mortensen AV: **Improved culture methods for isolation of Salmonella organisms from swine feces.** *Am J Vet Res* 2000, **61(11)**:1426-1429.
13. Marras SA, Kramer FR, Tyagi S: **Multiplex detection of single-nucleotide variations using molecular beacons.** *Genet Anal* 1999, **14(5-6)**:151-156.
14. Marras SA, Tyagi S, Kramer FR: **Real-time assays with molecular beacons and other fluorescent nucleic acid hybridization probes.** *Clin Chim Acta* 2006, **363(1-2)**:48-60.
15. Tyagi S, Kramer FR: **Molecular beacons: probes that fluoresce upon hybridization.** *Nat Biotechnol* 1996, **14(3)**:303-308.
16. Ellingson JL, Anderson JL, Carlson SA, Sharma VK: **Twelve hour real-time PCR technique for the sensitive and specific detection of Salmonella in raw and ready-to-eat meat products.** *Mol Cell Probes* 2004, **18(1)**:51-57.
17. Josefsen MH, Krause M, Hansen F, Hoorfar J: **Optimization of a 12-Hour TaqMan PCR-Based Method for Detection of Salmonella Bacteria in Meat.** *Appl Environ Microbiol* 2007, **73(9)**:3040-3048.
18. Malorny B, Bunge C, Helmuth R: **A real-time PCR for the detection of Salmonella Enteritidis in poultry meat and consumption eggs.** *J Microbiol Methods* 2007, **70**:245-251.
19. Malorny B, Made D, Teufel P, Berghof-Jager C, Huber I, Anderson A, Helmuth R: **Multicenter validation study of two blockcycler- and one capillary-based real-time PCR methods for the detection of Salmonella in milk powder.** *Int J Food Microbiol* 2007, **117(2)**:211-218.
20. Malorny B, Paccassoni E, Fach P, Bunge C, Martin A, Helmuth R: **Diagnostic real-time PCR for detection of Salmonella in food.** *Appl Environ Microbiol* 2004, **70(12)**:7046-7052.
21. Massi MN, Shirakawa T, Gotoh A, Bishnu A, Hatta M, Kawabata M: **Quantitative detection of Salmonella enterica serovar Typhi from blood of suspected typhoid fever patients by real-time PCR.** *Int J Med Microbiol* 2005, **295(2)**:117-120.
22. Moore MM, Feist MD: **Real-time PCR method for Salmonella spp. targeting the stn gene.** *J Appl Microbiol* 2007, **102(2)**:516-530.
23. Reynisson E, Josefsen MH, Krause M, Hoorfar J: **Evaluation of probe chemistries and platforms to improve the detection limit of real-time PCR.** *J Microbiol Methods* 2006, **66(2)**:206-216.
24. Shannon KE, Lee DY, Trevors JT, Beaudette LA: **Application of real-time quantitative PCR for the detection of selected bacterial pathogens during municipal wastewater treatment.** *Sci Total Environ* 2007, **382(1)**:121-129.
25. Chen W, Martinez G, Mulchandani A: **Molecular beacons: a real-time polymerase chain reaction assay for detecting Salmonella.** *Anal Biochem* 2000, **280(1)**:166-172.
26. Liming SH, Bhagwat AA: **Application of a molecular beacon-real-time PCR technology to detect Salmonella species contaminating fruits and vegetables.** *Int J Food Microbiol* 2004, **95(2)**:177-187.
27. Patel JR, Bhagwat AA, Sanglay GC, Solomon MB: **Rapid detection of Salmonella from hydrodynamic pressure-treated poultry using molecular beacon real-time PCR.** *Food Microbiol* 2006, **23(1)**:39-46.
28. Daum LT, Barnes WJ, McAvin JC, Neidert MS, Cooper LA, Huff WB, Gaul L, Riggins WS, Morris S, Salmen A, et al.: **Real-time PCR detection of salmonella in suspect foods from a gastroenteritis outbreak in Kerr county, Texas.** *J Clin Microbiol* 2002, **40(8)**:3050-3052.
29. Rahn K, De Grandis SA, Clarke RC, McEwen SA, Galan JE, Ginocchio C, Curtiss R 3rd, Gyles CL: **Amplification of an invA gene sequence of Salmonella typhimurium by polymerase chain reaction as a specific method of detection of Salmonella.** *Mol Cell Probes* 1992, **6(4)**:271-279.
30. Cortez AL, Carvalho AC, Ikuno AA, Burger KP, Vidal-Martins AM: **Identification of Salmonella spp. isolates from chicken abattoirs by multiplex-PCR.** *Res Vet Sci* 2006, **81(3)**:340-344.
31. Csordas AT, Barak JD, Delwiche MJ: **Comparison of primers for the detection of Salmonella enterica serovars using real-time PCR.** *Lett Appl Microbiol* 2004, **39(2)**:187-193.
32. Eyigor A, Carli KT, Unal CB: **Implementation of real-time PCR to tetrathionate broth enrichment step of Salmonella detection in poultry.** *Lett Appl Microbiol* 2002, **34(1)**:37-41.
33. Fey A, Eichler S, Flavier S, Christen R, Hofle MG, Guzman CA: **Establishment of a real-time PCR-based approach for accurate quantification of bacterial RNA targets in water, using Salmonella as a model organism.** *Appl Environ Microbiol* 2004, **70(6)**:3618-3623.
34. Fukushima H, Tsunomori Y, Seki R: **Duplex real-time SYBR green PCR assays for detection of 17 species of food- or water-borne pathogens in stools.** *J Clin Microbiol* 2003, **41(11)**:5134-5146.
35. Hein I, Flekna G, Krassnig M, Wagner M: **Real-time PCR for the detection of Salmonella spp. in food: An alternative approach to a conventional PCR system suggested by the FOOD-PCR project.** *J Microbiol Methods* 2006, **66(3)**:538-547.
36. Hoorfar J, Ahrens P, Radstrom P: **Automated 5' nuclease PCR assay for identification of Salmonella enterica.** *J Clin Microbiol* 2000, **38(9)**:3429-3435.
37. Khan AA, Nawaz MS, Khan SA, Cerniglia CE: **Detection of multidrug-resistant Salmonella typhimurium DT104 by multiplex polymerase chain reaction.** *FEMS Microbiol Lett* 2000, **182(2)**:355-360.
38. Malkawi HI, Gharaibeh R: **Multiplex PCR for the direct detection of Salmonella enterica from chicken, lamb and beef food products.** *J Basic Microbiol* 2003, **43(4)**:328-336.
39. Nam HM, Srinivasan V, Gillespie BE, Murinda SE, Oliver SP: **Application of SYBR green real-time PCR assay for specific detection of Salmonella spp. in dairy farm environmental samples.** *Int J Food Microbiol* 2005, **102(2)**:161-171.
40. Nucera DM, Maddox CW, Hoiem-Dalen P, Weigel RM: **Comparison of API 20E and invA PCR for identification of Salmonella enterica isolates from swine production units.** *J Clin Microbiol* 2006, **44(9)**:3388-3390.
41. Rychlik I, van Kesteren L, Cardova L, Svestkova A, Martinkova R, Sisak F: **Rapid detection of Salmonella in field samples by nested polymerase chain reaction.** *Lett Appl Microbiol* 1999, **29(4)**:269-272.
42. Wolffs PF, Glencross K, Norling B, Griffiths MW: **Simultaneous quantification of pathogenic Campylobacter and Salmonella in chicken rinse fluid by a flotation and real-time multiplex PCR procedure.** *Int J Food Microbiol* 2007, **117(1)**:50-54.
43. Wolffs PF, Glencross K, Thibaudeau R, Griffiths MW: **Direct quantitation and detection of salmonellae in biological samples without enrichment, using two-step filtration and real-time PCR.** *Appl Environ Microbiol* 2006, **72(6)**:3896-3900.
44. Kauffman F: **The diagnosis of Salmonella types.** Springfield Ill edition. 1950.
45. Altschul SF, Gish W, Miller W, Myers EV, Lipman DJ: **Basic local alignment search tool.** *J Mol Biol* 1990, **215(3)**:403-410.
46. Malorny B, Hoorfar J, Bunge C, Helmuth R: **Multicenter validation of the analytical accuracy of Salmonella PCR: towards an international standard.** *Appl Environ Microbiol* 2003, **69(1)**:290-296.
47. Lim YH, Hirose K, Izumiya H, Arakawa E, Takahashi H, Terajima J, Itoh K, Tamura K, Kim SI, Watanabe H: **Multiplex polymerase chain reaction assay for selective detection of Salmonella enterica serovar typhimurium.** *Jpn J Infect Dis* 2003, **56(4)**:151-155.
48. Soumet C, Ermel G, Rose N, Rose V, Drouin P, Salvat G, Colin P: **Evaluation of a multiplex PCR assay for simultaneous identification of Salmonella sp., Salmonella enteritidis and Salmonella typhimurium from environmental swabs of poultry houses.** *Lett Appl Microbiol* 1999, **28(2)**:113-117.
49. Soumet C, Ermel G, Rose V, Rose N, Drouin P, Salvat G, Colin P: **Identification by a multiplex PCR-based assay of Salmonella typhimurium and Salmonella enteritidis strains from environmental swabs of poultry houses.** *Lett Appl Microbiol* 1999, **29(1)**:1-6.
50. Carlson SA, Bolton LF, Briggs CE, Hurd HS, Sharma VK, Fedorka-Cray PJ, Jones BD: **Detection of multiresistant Salmonella typhimurium DT104 using multiplex and fluorogenic PCR.** *Mol Cell Probes* 1999, **13(3)**:213-222.
51. De Medici D, Croci L, Delibato E, Di Pasquale S, Filetici E, Toti L: **Evaluation of DNA extraction methods for use in combination with SYBR green I real-time PCR to detect Salmonella enterica serotype enteritidis in poultry.** *Appl Environ Microbiol* 2003, **69(6)**:3456-3461.
52. Herrera-Leon S, Ramiro R, Arroyo M, Diez R, Usera MA, Echeita MA: **Blind comparison of traditional serotyping with three multi-**

plex PCRs for the identification of Salmonella serotypes. *Res Microbiol* 2007, **158(2)**:122-127.

53. Pan TM, Liu YJ: **Identification of Salmonella enteritidis isolates by polymerase chain reaction and multiplex polymerase chain reaction.** *J Microbiol Immunol Infect* 2002, **35(3)**:147-151.
54. Pathmanathan SG, Cardona-Castro N, Sanchez-Jimenez MM, Correa-Ochoa MM, Puthucherry SD, Thong KL: **Simple and rapid detection of Salmonella strains by direct PCR amplification of the hilA gene.** *J Med Microbiol* 2003, **52(Pt 9)**:773-776.

Publish with **BioMed Central** and every scientist can read your work free of charge

"BioMed Central will be the most significant development for disseminating the results of biomedical research in our lifetime."

Sir Paul Nurse, Cancer Research UK

Your research papers will be:

- available free of charge to the entire biomedical community
- peer reviewed and published immediately upon acceptance
- cited in PubMed and archived on PubMed Central
- yours — you keep the copyright

Submit your manuscript here:
http://www.biomedcentral.com/info/publishing_adv.asp

